

CleanCast[™]

Advanced cast saw and dust extraction systems

introducing the Red Dot award winning design

QuietCast[™] paediatric saw system





CleanCast[™]

Oscillating cast saw systems

Contents

Cordless saw page 2

Paediatric saw page 3

Low voltage saws page 5

Mains voltage saw page 7

Accessories page 8

Saw blades page 9





De Soutter Medical CleanCast™ and QuietCast™ products comply with European Medical Device Directive 93/42 EEC and sets the highest standards for particulate filtration efficiency thus assisting conformance with COSHH and other international particle exposure regulations.

All products in this catalogue are supported by the De Soutter Medical fully comprehensive warranty.







Battery Powered Saw

Cordless

The De Soutter Medical CSB-100 saw allows the user total freedom of movement, which is particularly useful for medical emergency services, ward work, the veterinary profession and numerous other applications where the cordless concept is the only practical solution.

re the cordless concept is the rution.

CSB-100

Saw specification

Saw model	CSB-100			
Oscillating speed	16000 cpm			
Battery voltage	10.8 Vdc			
Weight	1050 grams			
Order code	1286224			

Supplied with saw blade Part No. 6250. See page 9 for other saw blade options.

Battery

Voltage	10.8 Vdc NiMH
Capacity	2600 mAh
Order code	15010



Battery charger

Supplied with interchangeable mains supply adaptors for use in most countries.

Charger model	BC-100
Supply voltage	100 - 240 Vac
Supply frequency	50 - 60 Hz
Order code	15020

CSB-100-A Cast Saw System

Order code 1286774 - comprising of:

CSB-100 Saw	Part No. 1286224
CB-100 Battery x 2	Part No. 15010
BC-100 Charger	Part No. 15020
Saw blade	Part No. 6250
Hexagon wrench	Part No. 204903
Carry case	Part No. 15030



Battery/Mains Powered Saw

Paediatric




Saw specification

Order code	1286994
Saw model	CSP-201
Oscillating speed	14800 - 17500 cpm
Handpiece weight	319 grams
Noise level	55 dB(A)

- Soft start feature with noise control helps minimise patient anxiety.
- · Rapid charging time under 3 hours.
- · Continuous 2 hour cutting performance.
- Built-in service reminder 300 charging cycles.
- Display allows the user to select the optimum cutting speed and view the remaining battery capacity.
- Micro processor controlled voltage booster ensures constant cutting performance - zero deterioration. See graph below.

CSP system power/cutting performance

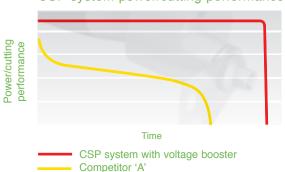

Battery specification

Order code	15250
Battery pack model	BP-201
Battery capacity	30.8 Wh
Battery weight	437 grams
Charging time from flat	< 3 hours

An integral clip with both battery pack or mains power adaptor for increased portability.





CSP-201-A Pro Cast Saw System

Order code 1286984 - comprising:

CSP-201 Saw and cable	1286994
BP-201 Battery Pack Pro	15250
IC-201 Inter Connector	15270
PS-201 Power Supply	15280
Carry Case	15290
Wrench	15300
Saw blade	15310
Saw blade	15320

Please specify Cordset required. See below

Mains Supply Cordsets

2 metres long

15610
15620
15630
15640
15650
15660

Low Voltage Saws

Corded electric

The CC5 and CC6 models combine the latest technological features with a functional, ergonomic handgrip. The smooth contours are designed to fit comfortably in the hand while providing the ideal cutting position. The balanced weight distribution increases manoeuvrability and helps to reduce wrist fatigue.

On the CC5 model, the dust extraction path is ducted internally through the handpiece which significantly improves visibility and handling in comparison to conventional saws fitted with extraction.

Rubber motor mounts insulate the internal mechanism from the outer case. This results in extremely low noise and reduces vibration to the operators hand.

A two speed control on the handpiece gives a quieter, low speed setting which is ideal for use with children and a high speed setting for fast removal of thick body casts.



Power Supply

CC5 and CC6 saws require their own low voltage power supply unit. On the CC5 model, the power supply is integral with the dust extraction unit. CC6 Saw is supplied by the PS6 power supply unit. See opposite page for options.



See opposite page for ordering instructions. See page 8 and 9 for accessories and saw blade options.

68 dB(A)

Sound level

Length

Power Supply for CC6

The power supply unit can be either bench mounted or wall mounted with the fixings provided.

Output voltage	22/29 Vdc
Supply frequency	50 - 60 Hz
Weight	1510 grams
Dimensions - mm	215L x 63H x 98W





The four stage filtration system consisting of an easy to change disposable collection bag, two cloth filters and a micro filter. The micro filter retains 99.997% of all dust particles greater than 0.5 microns.

Suction can be switched between powerful 'high' and quiet 'low' setting to suit operator and patient convenience.

The extractor with its exceptional dust filtration capability can be connected to either the CC5 or CC4 saw. When the saw is switched on, the extractor starts automatically. Likewise, the extractor stops automatically five seconds after the saw is switched off, thereby ensuring the system is purged of all cast particles.

The hose attachment position and low centre of gravity make the extractor very stable thus preventing accidental overturning.

Extractor Specification

Motor power (max.)	1 kW
Supply voltage	240 or 100 - 120 Vac
Supply frequency	50 - 60 Hz
Suction	2500 water gauge mm
Filtration	99.997% > 0.5 microns
Dust bag capacity	7.5 litres
Sound level	58 dB(A)

See page 8 for replacement filters and paper collection bag.

How to order CC5 system - with extraction

System component	S	U.K. plug 230 V	Euro plug 230 V	U.S. plug 120 V (CSA)	Aus plug 230 V
CC5 saw	Part No.	1270554	1270554	1280664	1270554
PS5 power supply	Part No.	609053	609053	609983	609053
Hose kit	Part No.	6400	6400	6400	6400
Extractor		1190	7020	7140	9490
Complete system ord	ering code	CC5A-230 UK	CC5A-230 EU	CC5A-120 US	CC5A-120 AUS

How to order CC6 system - less extraction

CC6 saw	Part No.	1270664	1270664	1280774	-	
PS6 power supply	Part No.	609073	609083	628703	-	
Complete system orde	ering code	CC6A-230 UK	CC6A-230 EU	CC6A-120 US	-	

CC5 and CC6 saws, power supplies and extractor are type tested to EN60601-1:1990, BS5724: Part 1: 1989, IEC 601-1: 1988. 120 V model is CSA approved.

Mains Voltage Saw

High power

The CC4 with its powerful 110 watt motor maintains rapid and efficient cutting on the thickest of body casts.

The knurled grip close to the saw blade provides the operator with excellent control and flexibility of movement.

There is an optional clear extraction hood which ensures good visibility when cutting. The hood is fully rotatable to suit the most convenient cutting position.

See extractor on page 6 for full extractor specification.

CC4 saw specification

Oscillating speed	14000 cpm (or 16000 cpm)
Voltage	230 Vac (or 120 Vac)
Weight	1300 grams
Sound level	70 dB(A)
Length	333 mm

See page 8 and 9 for accessories and saw blade options.



Clean Ca



How to order CC4 system - with extraction

System compone	ents	U.K. plug 230 V	Euro plug 230 V	U.S. plug 120 V	AUS plug 230 V
CC4 saw	Part No.	1273884	1273884	1273994	1273884
Hose kit	Part No.	78472	78472	78472	78472
Extractor	Part No.	1190	7020	7140	9490
Complete system	ordering code	CC4A-230UK	CC4A-230EU	CC4A-120US	CC4A-230AUS

How to order CC4 system - less extraction


		U.K. plug 230 V	Euro plug. 230 V	US plug. 120 V	AUS plug 230 V
CC4 saw	Ordering code	1280004	1270994	1280554	1281444

Accessories



Saw Blades

For CleanCast™ saws

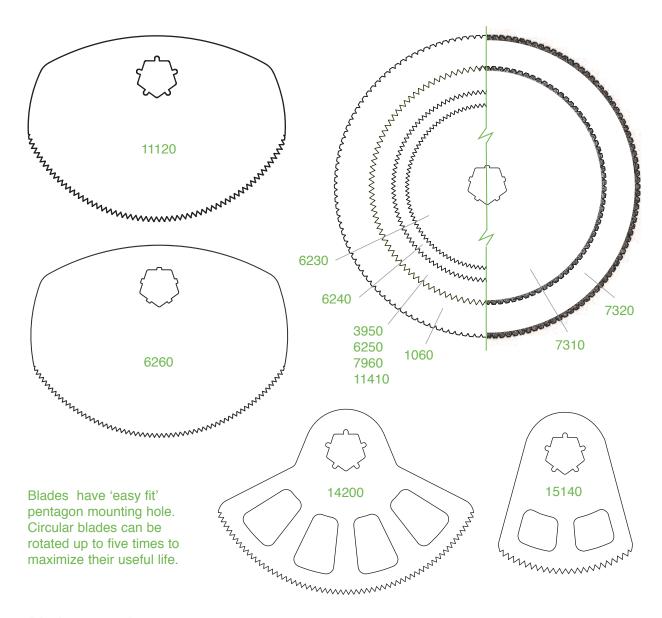


Blade range/Part no. All blades are sold in boxes of 5.

Blade dia mm	Blade shape	Hard chrome PTFE coated	Stainless steel	lon nitride	Titanium nitride	Diamond
44	Circular	6230				
50	Circular	6240				
64	Circular	6250	7960	11410	3950	7310
80	Circular	1060 1)				7320 1)
64	Segment		15140 ²⁾			
76	Segment		14200		11120	
80	Segment	6260				
POP - Gypsona		Optimal				
Fibreglass - Delt	alite/Scotchcast	Suitable				
Polyester - Delta	lite/Benicast/CellaCast	Not recommended				

¹⁾ Not suitable for Model CC5.

²⁾ Low noise/vibration blade - ideal for paediatric use.

Universal

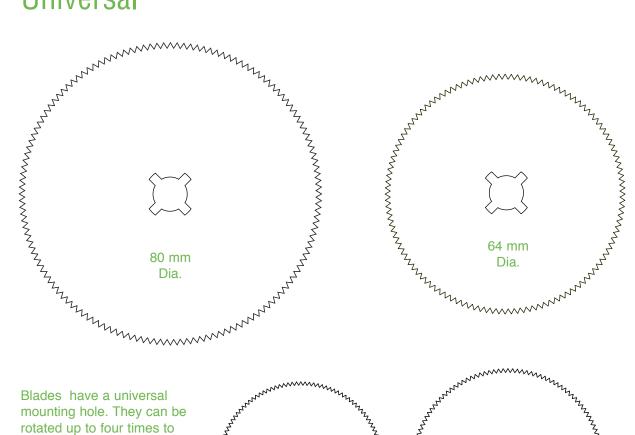



maximize their useful life. These blades can be used with these cast saw brands.

- Aesculap
- Zimmer
- M-Pact
- Martin
- BSN
- Medicon





Blade range/Part no. All blades are sold in boxes of 5.

Blade dia mm	Blade shape	Hard chrome PTFE coated	Titanium nitride
44	Circular	7040	
50	Circular	7800	
64	Circular	1180	7790
80	Circular	7810	
POP - Gypsona			
Fibreglass - Del	talite/Scotchcast		
Polyester - Delta	alite/Benicast/CellaCast		



For further information on De Soutter Medical cast removal systems and other De Soutter Medical products, please contact the address in your area

UNITED KINGDOM (INTERNATIONAL HEAD OFFICE)

De Soutter Medical Limited

Halton Brook Business Park, Weston Road Aston Clinton, Aylesbury, Bucks HP22 5WF

Tel: +44 (0) 1296 634 000 Fax: +44 (0) 1296 634 033 Email: info@de-soutter.com

AUSTRALIA

De Soutter Medical Australia Pty Ltd 2/12-14 Apollo Drive, Hallam

Victoria 3803

Tel: +61 (0) 3 9702 4441 Fax: +61 (0) 3 9702 4484

Email: australia@de-soutter.com

AUSTRIA

De Soutter Medical Austria Zweigniederlassung Österreich Dietrichsteingasse 10, A-3400 Klosterneuburg

Tel: +43 (0) 676 96 71 770 Fax: +43 (0) 2243 21 656 Email: austria@de-soutter.com

BELGIUM

De Soutter Medical Belgium In De Bruel 30 3620 Lanaken

Tel: +32 (0) 89/47 15 37 Fax: +32 (0) 89/70 12 19

Email: belgium@de-soutter.com

FRANCE

De Soutter Medical France 949 Avenue du Parc des Expositions 33260 La Teste de Buch

Tel: +33 (0) 5 56 54 89 36 Fax: +33 (0) 9 70 61 37 60 Email: france@de-soutter.com

www.de-soutter.com

GERMANY

De Soutter Medical Germany Bahnhofstraße 4 D-66625 Nohfelden

Tel: +49 (0) 68 52-99 12 46 Fax: +49 (0) 68 52-99 12 47

Email: deutschland@de-soutter.com

ITALY

De Soutter Medical Italia Località Fornace SNC 27022 Casorate Primo - PV Tel: +39 (0) 2 9009 4098

Fax: +39 (0) 2 9009 2673 Email: italy@de-soutter.com

UNITED STATES OF AMERICA De Soutter Medical USA, Inc. 224 Rolling Hill Road, Suite 12A Mooresville, NC 28117

Tel: +1 (704) 655 9040 Fax: +1 (704) 987 2035 Email: usa@de-soutter.com

NETHERLANDS

De Soutter Medical Nederland Gelderlandhaven 2X 3433 PG Nieuwegein

Tel: +31 (0) 85 0491480 Fax: +31 (0) 85 0491489

Email: nederland@de-soutter.com



113-FN-15