



hekaDrill




Ultimate Precision

The hekaDrill has been developed with innovation and precision as a guiding principle.

Key capabilities include:

- High speed motor control - up to 85,000 rpm. Allows high torque/low speed control. (min. 200 rpm)
- Advanced motor control ensures comfort for the surgeon during extended use with minimal heat gain.
- Adjustable bur exposure lengths.
- Wireless foot pedal operates using a secure encrypted network, ensuring greater user protection.
- Adjustable acceleration and braking functionality.

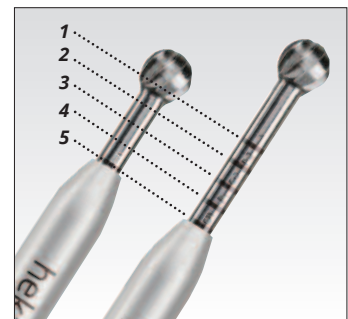

A single rotary bur allows 5 possible exposure positions offering total flexibility and reduced inventory costs

ELECTRIC HANDPIECE AND ATTACHMENTS

HANDPIECE	Model	Order code	Speed rpm
-----------	-------	------------	-----------

HDS-121
Cable length 4m

103347

200 - 85,000



BUR HEAD ATTACHMENTS	Model	Order code	Length - mm	Bur size - mm
STRAIGHT				



S-C07	103360	70	2.38
S-C08	103371	80	2.38
S-C10	103383	100	2.38
S-C15	103394	150	2.38
S-D09	103409	90	3.18
S-D14	103414	140	3.18

ANGLED



A-C07	103446	70	2.38
A-C10	103451	100	2.38
A-C15	103467	150	2.38
A-D14	103472	140	3.18

CRANIOTOME



C-C09 Small	103508	90	2.38
C-C10 Standard	103480	100	2.38
C-D12 Large	103495	120	3.18

PERFORATOR ATTACHMENT	Model	Order code	Speed rpm	Chuck type
-----------------------	-------	------------	-----------	------------

P-P01

103519

1,000

Hudson





Used with HDS-121 at 72,000 rpm

Special order only

See back page for range of cutters and burs



Description	Model	Order code
Console - less irrigation	PC-120	103305
Console - with irrigation	PC-121	103318

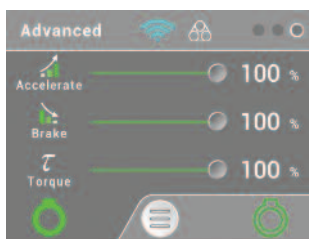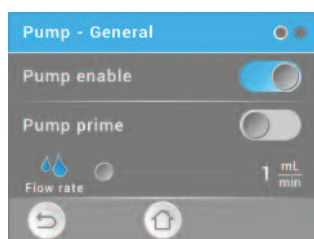

The advanced hekaDrill power console is available with or without irrigation and offers numerous features:

- Simultaneous dual handpiece operation.
- Adjustable acceleration and braking functionality.
- Integrated and customisable irrigation control with anti drip mode.
- 200+ user profiles can be stored and easily accessed.
- Rings illuminate to indicate port status for ease of use.

IRRIGATION AND FOOT PEDALS

STERILE IRRIGATION KIT



Sterile irrigation kit.
Box of 5.
Order code: S101461

FOOT PEDALS



- Both wired and wireless foot pedal available.
- Wireless foot pedal operates using a secure 128 bit encrypted network for complete protection.
- A or B power output selection.
- Foot pedal underlighting to allow for easy detection.



4 metre cable for
connecting and
charging wireless
foot pedal via
power console.

Order code: 103597

	Model	Order code
Wired	MFP-120	103556
Wireless	MFP-121	103562

POWER SUPPLY

Mains power cable 3 metres long. Please specify from below.

UK	Order code 8290
Euro	Order code 8300
USA/Canada	Order code 8310
Australia/New Zealand	Order code 8320
Japan	Order code 8330
Swiss	Order code 14700



HANDPIECE	Model	Order code	Speed rpm
-----------	-------	------------	-----------



HDS-120 103339 100,000

PERFORATOR	Model	Order code	Speed rpm	Operating pressure
------------	-------	------------	-----------	--------------------



P-P02 103524 850 7 bar

Uses Hudson chuck type

FOOT CONTROL	Model	Order code	Foot pedal inlet connection	Foot pedal outlet connection
	FP-120	103543	Integral 3 metre hose with MA7 probe	Accepts 'wash assist' hose from sizes below



3 metre length - Order code 103535
4 metre length - Order code 103803
5 metre length - Order code 103816



Wash cap for use with
HDS-120/ P-P02
Order code 103796

STERILISATION TRAYS AND CONTAINER

STERILISATION TRAYS AND CONTAINER



ELECTRIC INSTRUMENT SYSTEM

Order code 103610

Includes lid, container, electric tray and bur head tray.

Optional tray contents -

- 1 off HDS-121 handpiece
- 1 off P-P01 perforator
- 1 off electric cable

Bur head attachment tray

Optional tray contents -

- 5 off bur heads

PNEUMATIC INSTRUMENT SYSTEM

Order code 103629

Includes lid, container, pneumatic tray and bur head tray.

Optional tray contents -

- 1 off HDS-120 handpiece
- 1 off P-P02 perforator
- 1 off pneumatic hose

Bur head attachment tray

Optional tray contents -

- 5 off bur heads

Container and trays only

Wire container and lid

Tray for electric system

Tray for pneumatic system

Tray for bur heads

Order code 103944

Order code 103953

Order code 103989

Order code 103966

STERILE FILTERED CONTAINER

Accepts wire container plus wire trays.

Order code 103977

Lid gasket

Order code 669903

Filter paper - single use

Pack of 100 -

Order code 657243



Security tag

Pack of 100

Order code 657253




Sterile indicator label

Pack of 100

Order code 657263



ADJUSTABLE STERILE ROTARY BURS

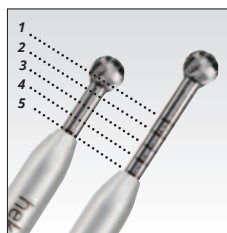

A single rotary bur allows 5 possible exposure positions offering total flexibility and reduced inventory costs



Fluted ball



Fine diamond ball



Coarse diamond ball



Extra coarse diamond ball



Fluted acorn



Fluted match



Fine diamond match



Extra coarse diamond cylinder



Fluted craniotome

Head size	Shaft dia.	Bur length - Order code						
		70 mm	80 mm	90 mm	100 mm	120 mm	140 mm	150 mm
2.0	2.38	C07-02021						C15-02021
3.0	2.38	C07-03021			C10-03021			C15-03021
4.0	2.38	C07-04021			C10-04021			C15-04021
5.0	2.38	C07-05021			C10-05021			C15-05021
6.0	2.38	C07-06021						C15-06021
3.0	3.18			D09-03021			D14-03021	
4.0	3.18			D09-04021			D14-04021	
5.0	3.18			D09-05021			D14-05021	
6.0	3.18			D09-06021			D14-06021	
1.0	2.38	C07-01051						
2.0	2.38	C07-02051						
4.0	2.38							C15-04051
5.0	2.38							C15-05051
4.0	3.18						D14-04051	
5.0	3.18						D14-05051	
6.0	2.38							C15-06071
6.0	3.18			D09-06071			D14-06071	
3.0	2.38	C07-03081						
4.0	2.38	C07-04081						
5.0	2.38	C07-05081						
5.0	2.38	C07-05023	C08-05023		C10-05023			C15-05023
6.0	2.38	C07-06023	C08-06023		C10-06023			C15-06023
7.5	2.38	C07-07523	C08-07523					C15-07523
6.0	3.18			D09-06023			D14-06023	
7.5	3.18			D09-07523			D14-07523	
9.0	3.18			D09-09023				
1.7	2.38	C07-01724	C08-01724		C10-01724			C15-01724
2.2	2.38	C07-02224	C08-02224		C10-02224			C15-02224
3.0	2.38	C07-03024	C08-03024		C10-03024			C15-03024
3.0	3.18			D09-03024			D14-03024	
3.0	2.38							C15-03054
3.0	3.18			D09-03054			D14-03054	
6.5	2.38							C15-06582
6.5	3.18			D09-06582			D14-06582	
1.5	2.38			C09-01535	C10-01535			
2.3	2.38			C09-02335	C10-02335			
3.0	3.18					D12-03035		

Special order only

Distributed by De Soutter Medical, manufactured by Zethon

213-EN-7

deSoutter
MEDICAL

De Soutter Medical Limited
Halton Brook Business Park, Weston Road, Aston Clinton
Aylesbury, Bucks HP22 5WF
Tel: +44 (0) 1296 634 000 Fax: +44 (0) 1296 634 033
Email: info@de-soutter.com
Web: www.de-soutter.com